
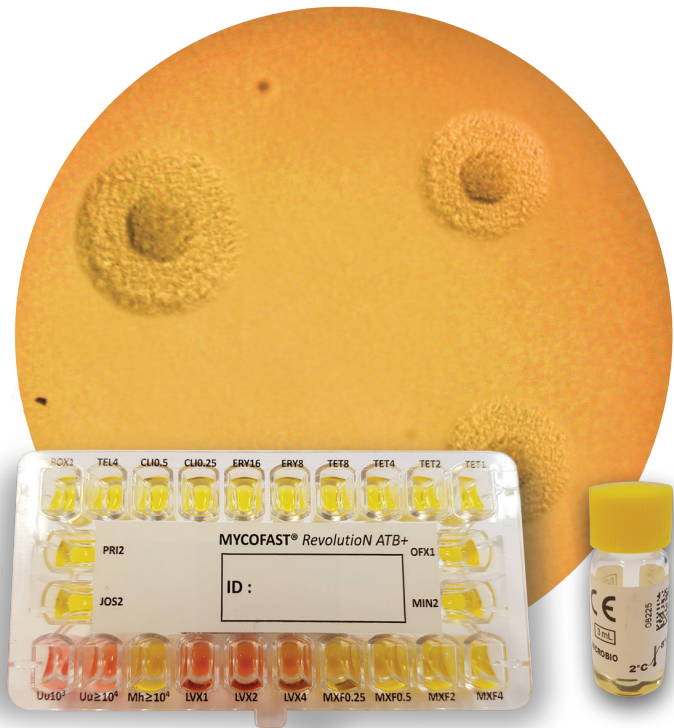


MYCOFAST[®] *Revolution* ATB+

UROGENITAL MYCOPLASMA DIAGNOSIS
DIRECT METHOD

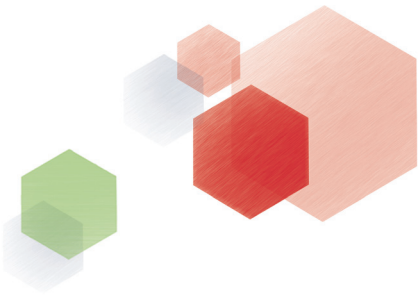

Expanding the scope of treatment options!

*Detection, Enumeration, Identification
and Antimicrobial Susceptibility Testing*

- ◆ Extended panel of antibiotics tested, in full compliance with worldwide recommendations
- ◆ Complete diagnosis in an all-in-1 product
- ◆ The right pathologic threshold tested for each specimen type
- ◆ Fast results for clinicians
- ◆ Visual, colorimetric reading

BENEFITS

Complete diagnosis in 24h*

- Detection, Enumeration
- Identification
- Antimicrobial Susceptibility Testing

Cutting-edge performance

Sensitivity

- Detection and enumeration of *U. urealyticum* at 10^3 and $\geq 10^4$ CCU/mL
- Detection and enumeration of *M. hominis* at 10^4 CCU/mL
- The right pathologic threshold tested for each specimen type

Specificity

- Identification secured by checking *U.u/M.h* natural resistance to antibiotics (Clindamycin and Erythromycin)

Reliability

- 6 antibiotics (TEL/LVX/MXF/ERY/TET/CLI) tested in compliance with CLSI guidelines (Clinical & Laboratory Standards Institute)
- 5 other antibiotics tested (MIN/OFX/ROX/PRI/JOS)

Fast results

- Integrated growth activator for *M. hominis*


Ease of use

- Ready-to-use reagents
- Each well pre-filled with growth medium
- Visual colorimetric reading
- No instrument required

Flexibility and safety

- Optimized transport medium: Specimens preserved up to 72 hours before analysis
- Clip-on closing system for the tray

METHODOLOGY

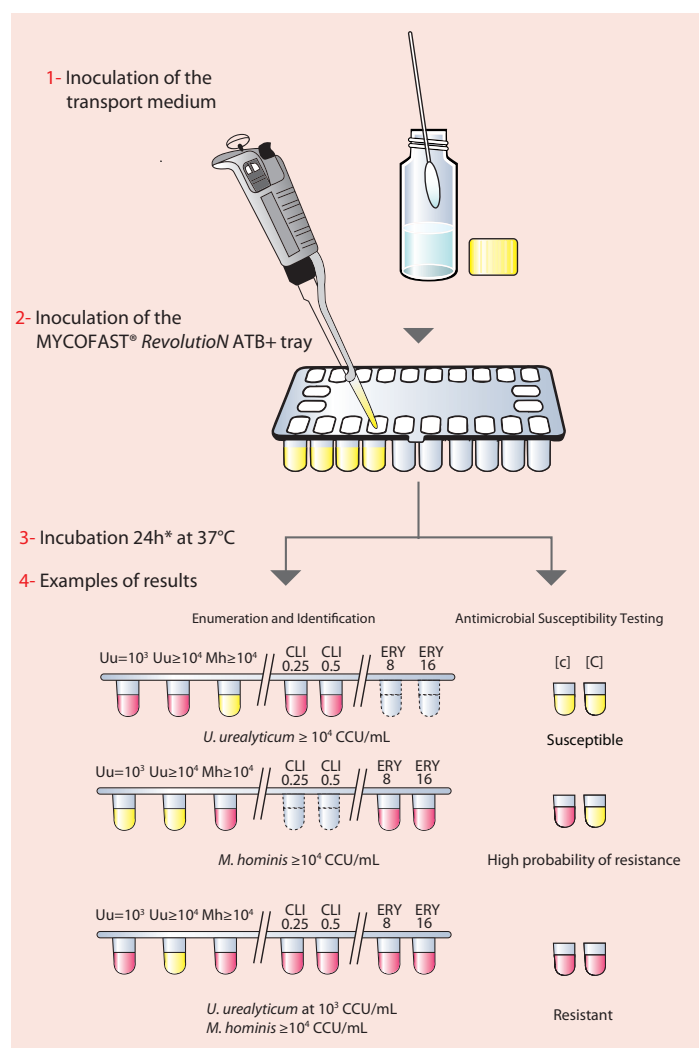

* Reading of all samples in 24h.
Second reading needed for liquid samples having negative results within 24h

PRODUCT RANGE

REF.	NAME	QTY	PURPOSE
00070	MYCOFAST® RevolutionN ATB+	25 tests	Complete kit
00061	UMMt Revolution (optional)	50 vials	Additional transport medium
00073	COMPLEMENT MYCOFAST® RevolutionN ATB+ (optional)	25 tests	Additional test trays
00064	MYCOPLASMA Stabilizer (optional)	100 tests	Additive for sample storage > 56 hours
00900	MYCOPLASMA Control (optional)	12 tests	U.u. and M.h. control strains

Product not available in the U.S
Please contact your sales representative for product availability in your country

 : ELITech Microbio (France)

WORLDWIDE OFFICES

Headquarters T: +33 1 41 45 07 10
Australia T: +61 1800 815 098
Belgium & Luxembourg T: +32 92 820 531
Brazil T: +55 27 3025 1415
France T: +33 4 83 36 10 82
Italy T: +39 02 48 40 35 42
New Zealand T: +64 0800 555 611

Serbia T: +381 11 2467119
Switzerland T: +41 26 663 86 60
The Netherlands T: +31 313 430 581
UK T: +44 1442 869320
United States T: +1 609 216 7361

 **ELITechGroup**
SOLUTIONS
tailored to your needs
www.elitechgroup.com
info@elitechgroup.com